



Free-standing room covering rail system

– Easy, practical and safe

- **Easy to install and remove**
- **Suitable for temporary lifting needs**
- **Suitable for reuse**
- **Height-adjustable – max. height 284 cm**
- **Requires no structural changes to the room or installation in walls or ceiling**
- **Does not take up any floor space and allows more room for the user and carer**

A mobile rail system that can be installed and removed quickly and easily as well as moved from place to place.

The system requires no structural changes to the room or installation in walls or ceiling, making it well-suited to temporary lifting needs in private homes, for example.

The free-standing rail system can easily be installed and removed by 1-2 people without the need for special tools.

The rail system is available in two standard dimensions, 3 x 3 metres and 4 x 4 metres, and is height-adjustable up to a maximum of 284 cm. Both GH1 and GH3 lifting modules can be used with the room covering system.



V. Guldmann A/S
Graham Bells Vej 21-23A
DK-8200 Aarhus N
Tel. +45 8741 3100
Fax +45 8741 3131
info@guldmann.com
www.guldmann.com

Free-standing room covering rail system, height adjustable

Technical specifications

Functions

Max. capacity : 350 kg

Weight and Materials

Free-standing room covering

rail system, 3 x 3 m : 151 kg
4 Vertical poles : 11 kg, each
4 B-rails : 16,4 kg, each
4 top corner brackets : 4,8 kg, each
1 Traverse rail
with trolleys : 22,2 kg

Free-standing room covering

rail system, 4 x 4 m : 178 kg
4 Vertical poles : 11 kg, each
4 B-rails : 21,8 kg, each
4 top corner brackets : 4,8 kg, each
1 Traverse rail
with trolleys : 27,6 kg

Rails and poles : Anodized aluminium
Brackets : Powder varnished steel

Measurements

Height min. : 2040 mm
Height max. : 2840 mm

Length, 3 x 3 m : 3200 mm
Width, 3 x 3 m : 3200 mm




Length, 4 x 4 m : 4200 mm
Width, 4 x 4 m : 4200 mm

Labelling

The product is manufactured in compliance with the Council Directive 93/42/EEC of June 14th 1993, including amendments, as medical device class 1.

Item Number/Model

555262 : Free-standing room covering rail system, height adjustable - 3 x 3 m
555570 : Free-standing room covering rail system, height adjustable - 4 x 4 m

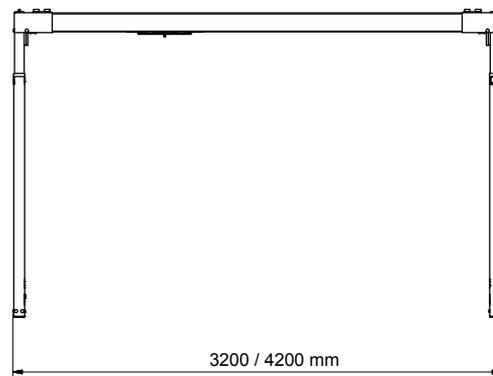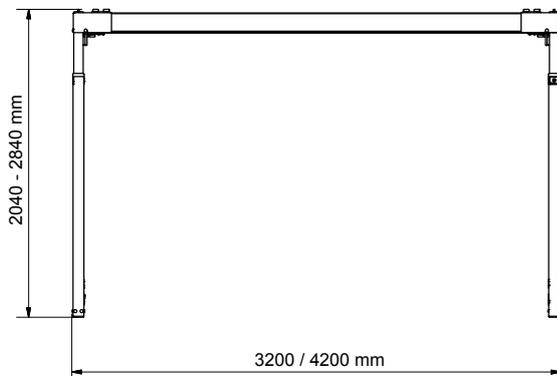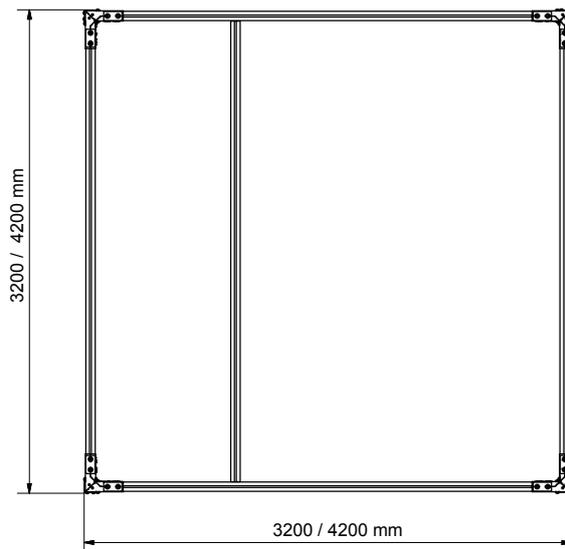

Time to care