

GH1

GH1 Q, GH1 F

Guldmann™

GH1 ceiling lift system



GH1 Q
– click on / click off version

Ceiling hoist for long term care

The GH1 ceiling lift system is designed to meet the full range of day-to-day lifting needs in homes and institutional environments. It takes up very little space, and is very easy to use with the full range of Guldmann lifting slings and other lifting accessories. The GH1 lifting module is also available in a GH1 Q version, which can easily be clicked-on and clicked-off the rail system. The GH1 and GH1 Q can be mounted in most types of rails.

High lifting quality and great comfort

GH1 combines the greatest possible comfort with efficient and safe working procedures. For the caregiver this means an improved working environment without needless manual lifts or strain – and for the user, a dignified lift with greater comfort and safety.

Easy to use

Just as the GH3 system the GH1 system enables the care staff to lift and move users without risking any back strain. The GH1 ceiling lift quickly and safely lifts users weighing as much as 205 kg (450 lbs), operated by an ergonomically correct and easy-to-use hand control. Quicker lifting frees up more time for highly important nursing and caring tasks.

Quick recharging

The GH1 is quickly recharged by placing the hand control in the GH1 recharging dock. Recharging normally only takes 2-3 hours, and a fully charged lifting module can cope with more than 55 lifts of 85 kg (187 lbs).

Click-on/click-off

The GH1 Q lifting module has a quick-release system that makes it easy to click on and off the rail to move it to another room. No tools are required!



GH1 ceiling lifts



Design that blends in

Both the GH1 ceiling lift and the rail system are easy to decorate so that they match or blend easily into any interior design of the space in which they are used. The GH1 is delivered with a choice of four different covers as standard.

More rails – more possibilities

The GH1 ceiling lift system is available with 3 different types of rails each specially designed for particular spans and lifting weights. The rails can easily be integrated with both ceiling and walls and can be customized to any interior design by coloured click-on covers. The covers can quick and easy be changed – which makes the rail system appear nice and newly painted – whenever a new resident moves into the room.

In addition to rails for mounting on wall and ceilings the GH1 system features a free standing rail system. These rails are easy to adjust in both height and width, providing an ideal way to deal with any temporary need for lifting help. The entire rail system can be set up by 1-2 people, with no tools needed. The entire system can be made ready for use – and taken down again – in only 10 minutes.

GH1 lifting modules can be mounted in most types of rails





GH1 ceiling lift sytem

Safe and reliable

The GH1 ceiling lift complies with all relevant standards, and is extremely safe and reliable in use.

The GH1 shares much of the same DNA as the award winning GH3 ceiling lift that has been installed in many of the most prestigious hospitals around the world. The GH1 ceiling hoist system makes full use of innovative and robust GH3 technology, design, functionality and costeffectiveness.

GH1 – ceiling lift

Userfriendly, powerful and fast lifting module for permanent installation and for installation in free-standing rail system.

- High lifting capacity – up to 205 kg (450 lbs)
- High lifting speed – 40 mm per sec with full load
- Operated by an intuitive and easy-to-use hand control
- Highly recyclable
- Extremely easy to clean
- GH1 Q – quick-release system makes it easy to relocate the lifting module from one room to another

GH1 – a ceiling lift that meets most demands

Free-standing rail system



GH1 F – the flexible ceiling lift



The GH1 F – is a flexible lifting module that is easy to move from one room to another by means of the transport and storage cart. It installs and removes itself automatically from the rail system, and eliminates any lifting, stretching or straining by the care staff in the process of moving and installing the lift.

With the GH1 F you can decide to install rails in several rooms and use fewer lifting modules that in the day-to-day care can easily be moved from one room to another.



Trolley 255 kg
max 255 kg
Guldmann

Guldmann

6001/07/15/16
max 255 kg
Guldmann

GH1 F – the flexible ceiling lift



Specially designed transport and storage trolley with space for the lifting module, the lifting hanger and the lifting sling.

Flexible, efficient and safe

Once installed the GH1 F gives you exactly the same functions as the GH1 lifting module – and delivers the same benefits as a fixed ceiling lift:

- Lifting height is not limited to the reach of the care staff
- User can be lifted from the floor or from chair to bed
- High lifting capacity – up to 205 kg (450 lbs)
- A wide range of detachable lifting hangers is available.

GH1 F – a flexible, self-installing ceiling lift

The GH1, GH1 Q and GH1 F are parts of the complete Guldmann range of lifting equipment that includes rails, ceiling lifts, lifting slings and accessories, – all designed as modules that work together perfectly. The modules can be combined in numerous different ways and can be installed just about everywhere.





The lifting module installs and removes automatically from the rail system – avoiding any strain to care staff.





US

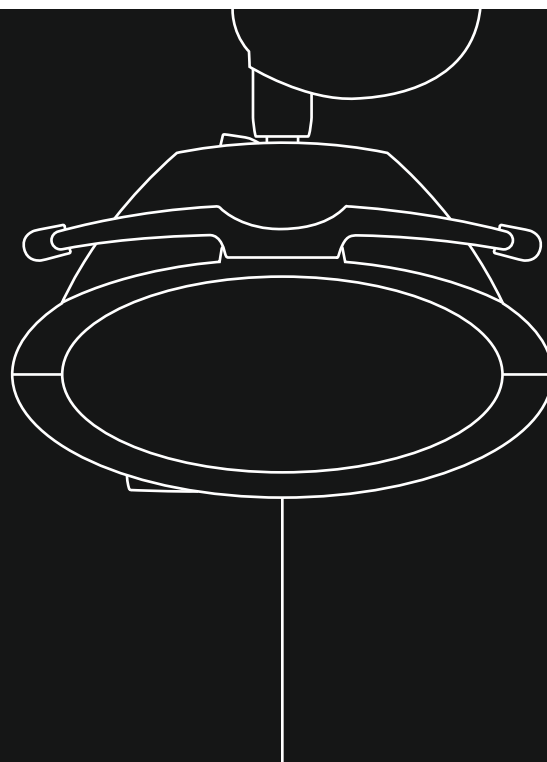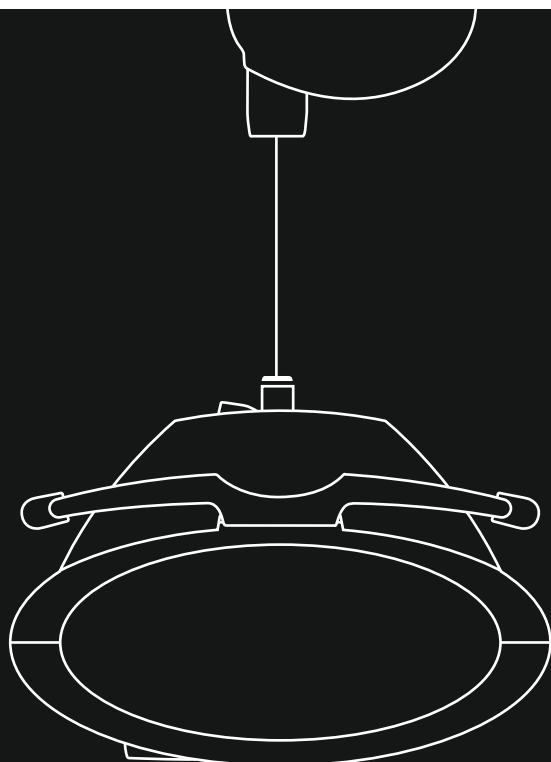

| Time to care |

V. Guldmann A/S
Corporate Office:
Graham Bells Vej 21-23A
DK-8200 Aarhus N
Tel. +45 8741 3100
info@guldmann.com
www.guldmann.com

Guldmann Inc.
14401 McCormick Dr.
Suite A
Tampa, FL 33626
Tel. +1 800 664 8834
Fax +1 813 880 9558
info@guldmann.net
www.guldmann.net

**Guldmann
Care-Lift Solutions ULC**
5716 Coopers Ave, Unit 9
L4Z 2E8 Mississauga, ON
Canada
Tel. +1 877 339 7779
Tel. +1 905 502 5438
ap@guldmann.net
www.guldmanncanada.ca

© Guldmann US-1830/03/2014

Guldmann™