



Wall Rail System

Freedom of movement – more time to care

***Discreet and easy to use and install
in both new and existing buildings***

The wall mounted rail system is a discreet system that creates a safe, comfortable and efficient solution for lifting and transferring of persons with loss of function. The system can also be used for gait training, balance training and for personal care tasks.

The wall rails are designed for direct installation on all kinds of walls and the flexible traverse rail will automatically compensate

for non-aligned walls of up to 80 mm (3 ¼ inches).

The wall-mounted rail system is used with either the GH3, GH1 / GH1 Q lifting module or the GHZ lifting module which is built into the traverse rail. When using the GH3, the lifting module is charged automatically in the rail, regardless of where the lifting module is placed in the rail system. The GH1, GH1 Q or GHZ are charged via the hand control in the charging station.

A continuous 'over the wall' rail system provides complete coverage and direct and efficient transfer from room to room e.g. from bedroom to bathroom. When not in use the system is almost invisible, as the lifting module, lifting hanger and sling can be stored in the cupboard. The lifting module can be reached from both rooms as the cupboard can be opened from both sides of the wall.

The wall rail is prepared for optional installation of LED light and can be used for hanging paintings, decorations, curtains and so on.



Discreet and elegant

- The system blends perfectly into the room and is almost invisible when not in use
- Integrated two-sided cupboard for storage of lift and slings
- Charging of lifting module either automatically in the rail (GH3) or via the hand control in a charging station (GH1/ GH1 Q, GHZ)
- The wall rails are white (RAL 9016) and blend discreetly into most environments




Simple and flexible system

- 'Over the wall' rail system
- Complete coverage and effective transfer from room to room
- Consecutively, comfortable and smooth movement
- Easy to operate



Easy and efficient

- Easy and economical to install
- Easy to remove and re-use
- Installation directly on all kinds of walls
- Flexible traverse rail automatically compensate for non-aligned (oblique) walls with up to 80 mm (3 1/4 inches)

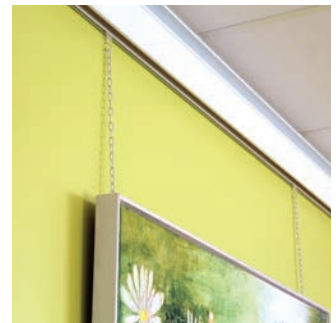

Extra functions

- Prepared for optional installation of LED lights
- Mounting of paintings, decorations, curtains, etc.

LED lights, suspension cords and hooks are not available in the Guldmann product range.

| Time to care |

Guldmann Inc.
USA
Tel. 800 664 8834
Tel. 813 880 0619
info@guldmann.net
www.guldmann.net

**Guldmann Care-Lift
Solutions ULC**
Canada
Tel. +1 877 339 7779
Tel. +1 905 502 5438
ap@guldmann.net