



GTP1

Transfer platform

Easy, comfortable and safe transfer in a standing position

The Guldmann GTP1 Transfer platform is a manual mobile working tool that facilitates moving standing and semi-independent users.

By actively participating in the transfer, the user is able to utilise and train his/her own muscle function and balance.

The platform is robust and stable, ensuring support and safety during short standing transfers. If needed, a transfer belt can be used for additional support during the standing transfer.

Typical moving and handling tasks using the Transfer platform include:

- To and from the toilet/commode chair
- To and from a shower chair
- Practising getting up/sitting down (mobilisation)
- Maintaining a standing posture (mobilisation)

Easy and safe to manoeuvre

The Transfer platform is easy to move around and requires very little space. It can be pushed, pulled and rotated around its own axis. This makes it very easy to manoeuvre, even with limited room.

Flexible and efficient


The Guldmann Transfer platform can be used in all healthcare environments and ensures an active, safe standing transfer of the user and a safe and ergonomic working environment for the carer.

The platform is easy to transport between different locations, for example in home care, as the top and bottom sections can easily be separated and transported in two parts.

Max. weight capacity of 175 kg.



Product details

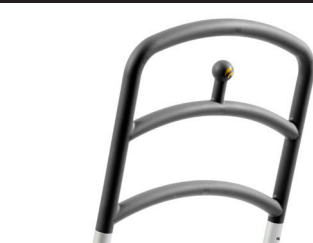

Multi-grip handle

The handle, which has a soft and non-slip surface, is ergonomically designed with curved grips on three levels. This allows both the user and the carer to choose the best place to hold onto, and the three handles can be used as a 'ladder' when the user stands up or sits down.



Low entry height

The baseplate is made of moulded aluminium with a non-slip surface. The entry height of just 30 mm and the shape ensures that it can be positioned very close to a wheelchair or toilet, for example. The low entry height makes it easy for the user to get their feet onto the baseplate.



Adjustable lower leg support

The adjustable lower leg support has a soft padding that provides good and comfortable support for the user. To achieve maximum support, the lower leg support can be easily and quickly adjusted in height and rotated (tilted) using the easy-grip function.



Central brake

The foot brake is centrally positioned and locks both wheels simultaneously. The brake is colour-coded to indicate clearly whether it is locked or not.



Transfer belt

For extra support and safety during transfer. It is placed around the user's upper body and attached to the central grip on the Transfer platform's multi-grip handle. NOTE: the Transfer belt should not be used to rise the user into a standing position.



Prepared for recycling

Smooth surfaces and solid construction make the GTP1 Transfer platform easy to clean and suitable for recycling. The Transfer platform is prepared for a simple 'end of life' separation and recycling of materials.

Time to care

V. Guldmann A/S
Tel. +45 8741 3100
info@guldmann.com
www.guldmann.com